

Fabric Thickness Measurement under Tensile Forces

Diana Šimić Penava^{1, a*}, Željko Penava^{2, b} and Željko Knezić^{3, c}

¹Faculty of Civil Engineering, University of Zagreb, Kačićeva 26, 10000 Zagreb, Croatia

^{2, 3}Faculty of Textile Technology, University of Zagreb, Pr. b. Filipovića 28a, 10000 Zagreb, Croatia

^adianas@grad.hr, ^bzpenava@tff.hr, ^czknezic@tff.hr

Keywords: thickness, woven fabrics, tensile force, extension.

Abstract. Woven fabric thickness (t) is its third dimension, which is considerably smaller than the fabric width and length, but at the same time it is very important because it affects the usage and properties of the fabric. During stretching, the fabric thickness is changing depending on the value of the tensile force. In an effort to be as precisely as possible to measure dynamic thickness changes of the fabric, and simultaneously collect all the required values about the change of the tensile force, extension and the area on which tensile force acts, new measuring apparatus are designed and manufactured. This measuring apparatus allows the simultaneous measurement of dynamic thickness of fabric, the fabric contraction and related tensile axial forces.

Introduction

Under fabric thickness the greatest distance between the surface of the front side and the mutually parallel surface of the back side of the fabric is meant. Fabric thickness is measured in mm. Fabric thickness plays an important role in manual processing, design and end-use of textiles. Fabric thickness change, depending on the effect of various parameters [1]. In order to clarify the effects of some specific parameters on fabric thickness, such as yarn count change, change of warp and weft density, change of warp and weft thread tension in the weaving process, it is necessary to carry out experimental tests of fabric thickness change. To allow a real understanding of thickness behavior when a fabric is stretched, it has been necessary to develop a precise new device for measuring fabric thickness. A textile fabric is a very flexible material and it is easy to imagine that, during such a test, the sample moves a lot in the three spatial directions. Moreover, a textile fabric is made of interlaced yarns constituted of fibers twisted together [2]. Various apparatuses for measuring static thickness have been developed. Haven used a new fabric thickness-measurer [3], Shiefer used a compressometer for thickness and resilience measurements [4]. The non-contacting measurements [5, 6] (optical) were too heavy to be fixed on the fabric without deforming it. It came out that the fabric should carry the sensor during the test.

The aim of this work was to measure the dynamic change of fabric thickness using a newly constructed device when the fabric was in the axial state of stress in the weft and warp direction.

Experiment

The experimental part deals with measurements of fabric thickness if tensile force acts on fabric samples until rupture. For this purpose a new device for dynamic measurements of

N is used. The electric motor work (1) is computer-controlled using the power and control unit with manual settings (12). In the middle of the test specimen a digital thickness gauge (8) with displacer (7) resting on the fabric specimen is placed. The displacer diameter is 2.54 mm, and the displacer area is 5.07 cm². The accuracy of thickness measurement is 0.001 mm. The full measurement range of the thickness gauge ranges from 0 to 25.4 mm. The displacer rests on the fabric under which the reference axle of the thickness gauge (5) is affixed to the support of the thickness measurement system (6). The applied digital thickness gauge is also furnished with displacement force setting (11). The whole system is connected to the computer (14) through the computer interface with A/D converter (13) and controlled by the serial port.

Results and discussion

For each fabric (F1 and F2), results are shown on diagrams Figures 2-3. The measures have been made on 2 samples in each direction in warp and weft. The Microsoft Excel software was used for statistical analysis of data at p<0.05 for five measurements.

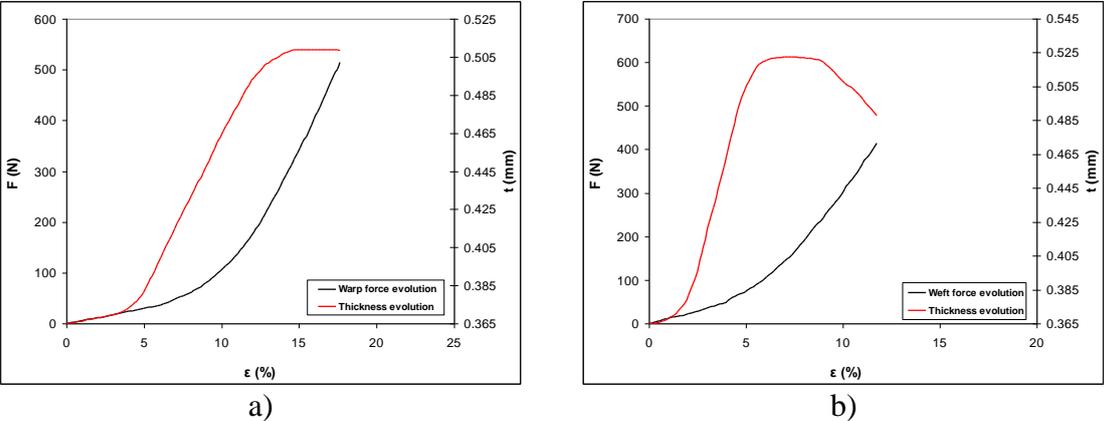


Fig. 2 Diagram of tensile force- relative extension – thickness (F- ε- t) for fabric F1: a) force acts in the warp direction, b) force acts in the weft direction

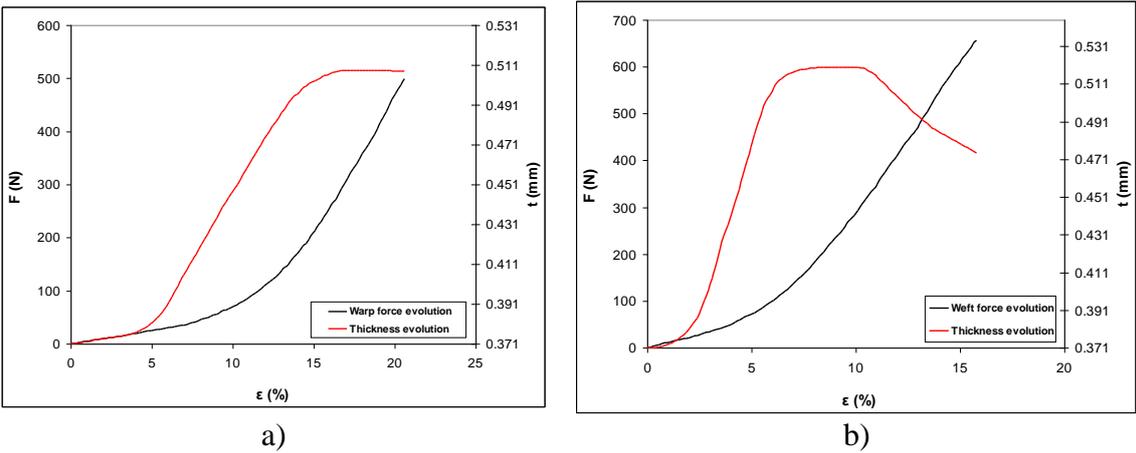


Fig. 3 Diagram of tensile force- relative extension – thickness (F- ε- t) for fabric F2: a) force acts in the warp direction, b) force acts in the weft direction

The known mechanisms of the fabric longitudinal deformation during a tensile test can be described as follows: first, low decrimping of yarns, and fabric structure rearrangement due to internal stresses; these two mechanisms occur under low loads and provide low decrimping of

stretched yarns; weave reaction and high de-crimping of the yarn in the direction of the test; and direct yarn strengthening when decrimping is complete. Before analyzing each fabric behavior, it is worth noticing the different initial thickness values, Table 1. These differences were expected because of the difference in fabric weft density. Fabric thickness measurement is highly dependant on the applied tensile force during the measure. Fabric thickness changes while the tensile force stretches the fabric. Thickness - axial force -extension curves in Figures 2 and 3 show typical non-linear relationships for fabric samples if tensile force acts in the weft or warp direction.

As is shown in Table 1, the number of yarns/cm is higher in warp than in weft. The woven fabric F1 has been tested because it is lightweight. Between 0 and 4% of strain in warp direction, and 0 and 2% of strain in weft direction, increasing thickness is very small. Between 4 and 12% of strain thickness is growing up from 0.375 to 0.505 mm within the warp. From 2% to 5% of strain in weft direction, thickness increases very fast from 0.370 to 0.525 mm.

For woven fabric F2, if the force acts in the warp direction in the extension range from $\epsilon=0\%$ to $\epsilon=4\%$ fabric thickness slightly increases to 0.381 mm. While the values of relative extension rise from $\epsilon=4\%$ to $\epsilon=15\%$, fabric thickness rapidly grows and reaches 0.511 mm. From 2% to 6% of strain in weft direction, thickness increases very fast from 0.379 to 0.529 mm.

Figures 2b and 3b show that in case if force acts in the weft direction with increasing weft density, the values of tensile forces and associated elongations at break of the fabric test specimen increase.

Conclusion

Measurements of fabric dynamic thickness were performed on a newly constructed measuring apparatus which simultaneously records all necessary values of changing tensile force and associated extension. Fabrics of different weft densities reach approximately the same maximum values of relative fabric thickness with almost equal forces and relative extensions. It follows that weft density does not affect fabric thickness if it is subjected to the action of tensile force in the weft direction. Test results have shown that an increase in axial tensile forces on the fabric causes an increase in fabric thickness for all fabric densities.

References

- [1] B. Olofsson, A general model of a fabric as a geometric-mechanical structure, *Journal of the Textile Institute Transactions*. 55 (1964) 541-557.
- [2] X. Xiao et al., Through-thickness permeability modelling of woven fabric under out-of-plane deformation, *Journal of Materials Science*. 49 (2014) 7563-7574.
- [3] G. B. Haven, New Fabric Thickness-Measurer, *Textile Research Journal*. 3 (1933) 229-237.
- [4] H. F. Schiefer, The compressometer: An instrument for evaluating the thickness, compressibility and compressional resilience of textiles and similar materials, *Textile Research Journal*. 10 (1933) 505-513.
- [5] A. Majumdar, S. S. Saha, A novel fabric thickness and compression tester using magnetic inductance principle, *Journal of the Textile Institute*. 99 (2008) 359-362.
- [6] J. Fan et al., A portable tester for nondestructively measuring fabric properties, *Textile research Journal*. 72 (2002) 21-26.